


LabSen 761 Electrode

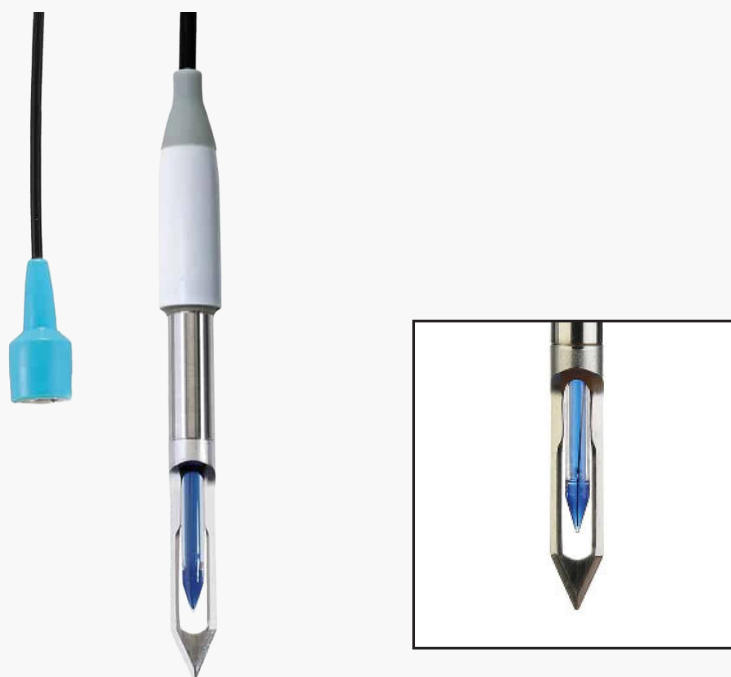

Main Features

- Integrated temperature sensor, patented Swiss sensor technology and materials
- Blade for pH measurement of meat, fish, and semi-frozen foods
- Food-grade titanium alloy blade, easy to clean
- Pointed glass membrane delivers stable readings in 15–30 seconds, up to 10 times more robust than conventional electrodes
- Polymer electrolyte with pinhole junction structure prevents clogging; no refilling required
- Long-life reference system (glass tube, AgCl, silver reference wire) improves stability and service life
- Integrated temperature sensor enables automatic temperature compensation
- Compatible with all pH meters with a BNC connector

Technical Parameters

Measuring range	(0 ~ 14) pH	Electrolyte	Polymer
Recommended operating temperature	(0 ~ 80)°C	Storage solution	3M KCl
Shaft material	Titanium blade	Membrane impedance	<250 MΩ
Glass membrane type	S	Electrode dimensions	(Ø12×115) mm
Reference system	Long-life	Connector	BNC
Junction	Ceramic + single pore	Cable	Ø3 mm×1 m